
waterpik®

Waterpik® Sensonic® Professional Toothbrush

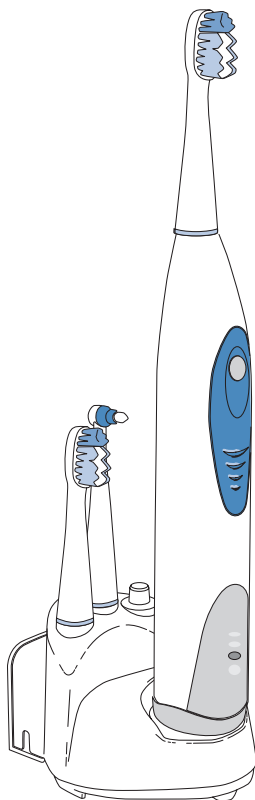**Waterpik® Sensonic®
Professional Zahnbürste**

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.....
www.waterpik.com

WELCOME

Congratulations, as an owner of the Waterpik® Sensonic® Professional Toothbrush, you are now on your way to better plaque removal, naturally whiter teeth, and a healthier smile!

The Sensonic® Professional Toothbrush utilizes an ultra high bristle speed and wide range of brushing motion. As a result, the Sensonic® Professional Toothbrush generates superior cleaning action and delivers improved plaque removal, healthier gums, and a fresh clean feeling.

Thank you for choosing the Waterpik® Sensonic® Professional Toothbrush.

TABLE OF CONTENTS

| | |
|--|----|
| Important Safeguards | 4 |
| Product Description | 6 |
| Getting Started | 8 |
| Care of your Waterpik® Sensonic® Professional Toothbrush | 10 |
| Limited Two-Year Warranty | 11 |

IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

As with any electrical product, especially when children are present, always follow these basic safety precautions. The safeguards below should always be followed.

PLEASE READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Do not handle the plug with wet hands.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not immerse charger in water or other liquid.
- Do not handle charger base while bathing.
- Do not reach for a charger that has fallen into any liquid. Unplug immediately.

WARNING:


To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:

- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it has dropped into any liquid. Return the product to a Water Pik, Inc. Authorized Service Center.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Be sure only the bristles of the toothbrush touch your teeth and gums when using the device.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- For household use only. Do not use outdoors.
- Use only those attachments, including toothbrushes and interdental brushes, which are recommended by Water Pik, Inc.
- Keep charging unit and cord away from heated surfaces.
- Do not operate or use where oxygen or aerosol sprays are being used.


- To ensure that your charger base does not accidentally fall and submerge in water, mount it on the wall.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult with your dentist before using this instrument or any other oral hygiene aid.
- Pacemaker wearers and persons with medical implants containing magnets or electronics: contact your physician before using this product.

SAVE THESE INSTRUCTIONS

PRODUCT DESCRIPTION

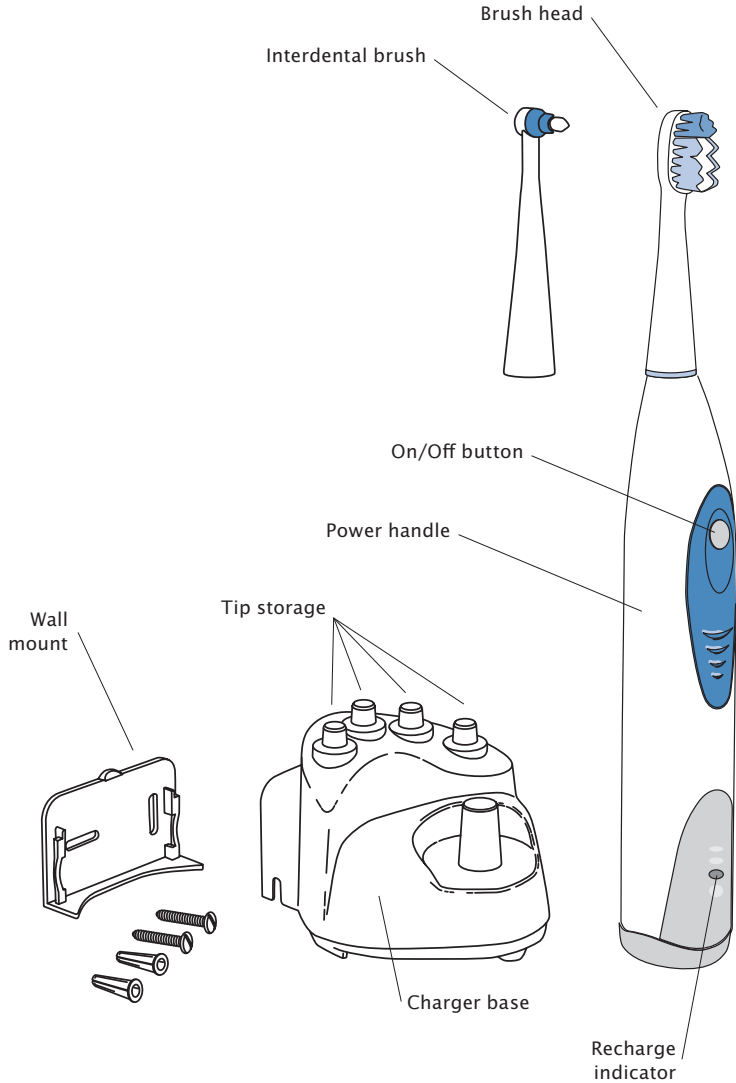
Your Sensonic® Professional
Toothbrush can be easily mounted
to the wall using the enclosed items.

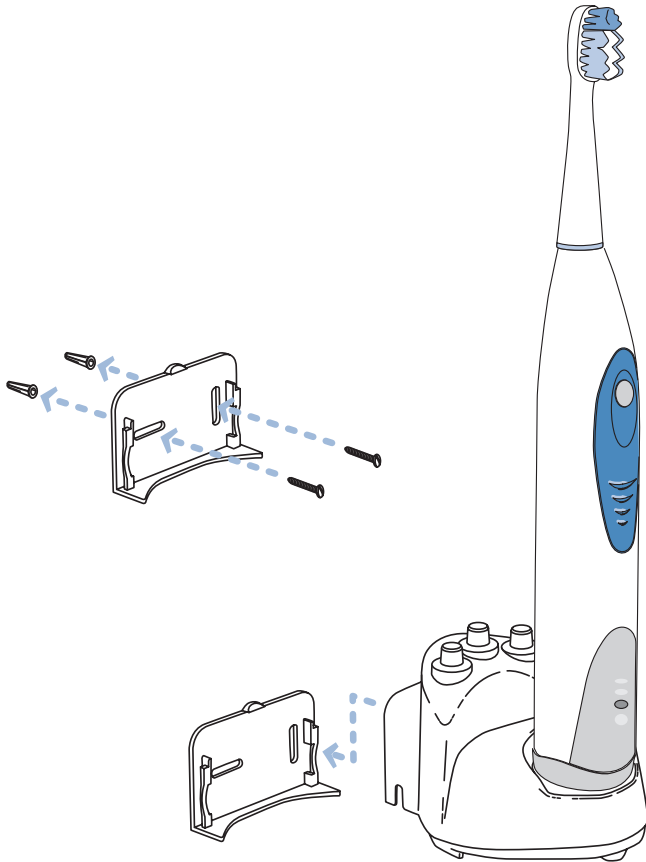GETTING STARTED


Charging the Sensonic® Professional Toothbrush

Charge for 16 hours before first use. The charge indicator on the handle will flash when charging and will stay on continuously when the unit is fully charged. A full charge will last for approximately 1 week of typical use (2 minutes of brushing twice a day).

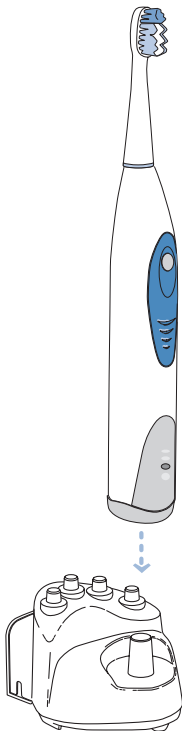
Sonic brushing

- Slide the brush head onto the handle. Push brush head down until it locks into place on handle.
- Place the brush head into your mouth before turning the unit on.
- Select low or high brushing speed by pressing power button after unit is on. Press button once for low speed, and twice for high speed.
 - Low speed – gentle cleaning and gum massage
 - High speed – superior plaque removal
- For best results, position brush along gum line at a 45-degree angle.
- Use a lighter brushing pressure than when using a manual toothbrush to optimize sonic cleaning effectiveness.

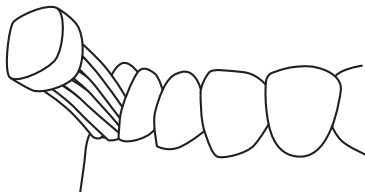- Move the brush slowly around entire mouth. Many dentists recommend spending 30 seconds in each of the four quadrants of your mouth.
- Quad Timer: your Sensonic® Professional Toothbrush will create a short stuttering sound every 30 seconds, to remind you to move to the next section of your mouth.

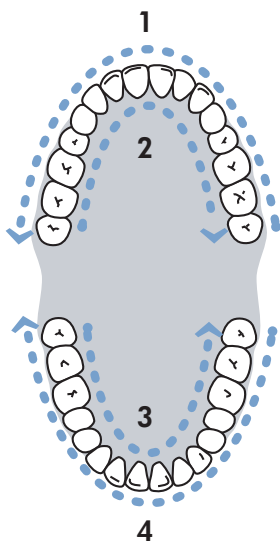Using the interdental brush

- Slide the interdental brush head onto the handle. Push brush head down until it locks into place.
- Brush may be used with or without toothpaste.
- Place brush tip perpendicular to gum line between two teeth.
- Turn on power handle and lightly follow crevices between teeth.

CARE OF YOUR WATERPIK® SENSONIC® PROFESSIONAL TOOTHBRUSH

1. After each use, turn off unit, remove the brush head from the handle and rinse brush head under water.
2. Wipe the handle dry after each use.
3. Periodically, clean the charger base by unplugging the unit and wiping with a damp cloth.

Service

There is no required routine service maintenance and this product cannot be repaired.

For all your warranty and accessory needs, please visit www.waterpik.com, click on "Product Support" and the "Service" link for your nearest Authorized Service Center. Shipment of product directly to Waterpik offices may delay service.

When outside of the United States: Contact your local dealer or write Water Pik, Inc. for the address of the nearest authorized service center.

Battery Disposal

The power handle contains a Nickel Metal Hydride battery, which does not contain heavy metal and is environmentally friendly. The battery is not replaceable. At the end of the battery life, please recycle the power handle at your nearest recycling center.

Important Precautions

Using Arm & Hammer® PeroxiCare® toothpaste, abrasive or chemical cleaners may damage the plastic housing of this product. Use of these products voids this product warranty.

LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will repair or replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as brush heads/ user replaceable batteries/etc.

This warranty gives you specific legal rights and you may also have other rights which vary from country to country.

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©2006 Water Pik, Inc.
Arm & Hammer® PeroxiCare® are registered trademarks of Church & Dwight Co, Inc.
Waterpik®, Waterpik® (stylized) and Sensonic® are registered trademarks of Water Pik, Inc.
Arm & Hammer® PeroxiCare® sind eingetragene Waren der Church & Dwight Co, Inc.
PeroxiCare® Arm & Hammer® sont des marques déposées de Church & Dwight Co, Inc.
Arm & Hammer® PeroxiCare® zijn gedeponeerde handelsmerken van Church & Dwight Co, Inc.